

Europäisches Patentamt  
European Patent Office  
Office européen des brevets



(11) **EP 1 048 950 A1**

(12) **EUROPEAN PATENT APPLICATION**

(43) Date of publication:  
02.11.2000 Bulletin 2000/44

(51) Int Cl.7: **G01N 33/566, G01N 30/46**

(21) Application number: **99201301.1**

(22) Date of filing: **26.04.1999**

(84) Designated Contracting States:  
**AT BE CH CY DE DK ES FI FR GB GR IE IT LI LU  
MC NL PT SE**  
Designated Extension States:  
**AL LT LV MK RO SI**

(72) Inventors:  
• van der Greef, Jan  
3971 BM Driebergen (NL)  
• Irth, Hubertus  
1078 NP Amsterdam (NL)

(71) Applicant: **ScreenTec B.V.**  
**2333 CC Leiden (NL)**

(74) Representative: **Ottevangers, Sietse Ulbe et al**  
**Vereenigde,**  
**Postbus 87930**  
**2508 DH Den Haag (NL)**

(54) **Identification of ligands for orphan receptors using on-line coupling of mass spectrometry to continuous-flow separation techniques**

(57) The present invention relates to the on-line coupling of mass spectrometry (MS) to continuous-flow separation techniques for detecting 'orphan' analytes, viz. receptors which do not have known ligands. In a fur-

ther embodiment, this on-line detection method is used as either a screening method for assays with an unknown ligand. Furthermore, this invention relates to compounds detected by this method and the use of these compounds as a ligand for affinity molecules.

Best Available Copy

EP 1 048 950 A1

## Description

[0001] The present invention relates to the on-line coupling of mass spectrometry (MS) to continuous-flow separation techniques for detecting 'orphan' analytes, viz. affinity molecules which do not have known ligands, such as orphan affinity proteins, e.g. orphan receptors. In a further embodiment, this on-line detection method is used as either a screening method for assays with an unknown ligand. Furthermore, this invention relates to compounds detected by this method and the use of these compounds as a ligand for affinity molecules.

[0002] A range of biochemical assays based on continuous flow biochemical detection has been developed during the recent years. Typically, the affinity interaction between ligand and affinity protein has been monitored by measuring the concentration of highly sensitive, usually fluorescence labeled reporter molecules using either heterogeneous or homogeneous biochemical detection strategies. When dealing with 'orphan' receptors however, i.e. receptors for which 'affinity' ligands have not been found so far, new strategies have to be developed, which enable so-called 'label-free detection' of affinity ligands. In contrast to the label-based assays, which usually allow the non-binding analytes from the sample to co-elute with the reporter molecules, label-free assays for orphan receptors demand strict separation of both binding and non-binding ligands.

[0003] By configuring restricted access columns and hollow fiber modules in the biochemical part of the system, affinity ligands can be isolated from complex biological matrices. Subsequently, coupling the biochemical setup to mass spectrometry allows the on-line detection of ligands exhibiting affinity for the orphan affinity molecule.

[0004] It is known from the prior art that biochemical assays are highly sensitive detection techniques, which combine the selectivity of biospecific interactions with the sensitive detection of labels used as reporter molecules.

[0005] An example of such a technique is described in WO 91/13354. This document describes a flow immunosensor wherein an antibody specific to an affinity molecule is immobilized on a support. The antibody's binding sites are saturated with a labeled form of a known ligand. A liquid to be analyzed is contacted with the antibody to allow any affinity molecule present in the liquid to displace the labelled antigen. Finally, the displaced labelled antigen is detected.

[0006] However, it is known as well that biochemical assays suffer from the problem of cross-reactivity, which means that antibodies react with more than one analyte. This leads to erroneous results. For this reason, biochemical assays are frequently combined with a fractionation step, e.g., a separation step using HPLC or another type of liquid chromatography.

[0007] In the field to which the present invention relates, there is a demand for on-line coupling of a frac-

tionation step and an biochemical assay detection system. Several approaches have been proposed and described to perform continuous-flow biochemical assays. Most of these biochemical assays are in the form of a postcolumn reaction detection system based on a sequential addition type biochemical assay.

[0008] Cassidy et al., *Anal. Chem.* **64** (1992), 1973-1977 disclose a kinetically controlled immunoassay based on the sequential addition of antibody, sample and label on a protein A column, performing immunoassays for albumin and transferrin in less than 1 minute.

[0009] Nilson et al., *J. Chromatography*, **597** (1992) 383-389 describe a continuous-flow competitive assay involving enzyme-labelled antibodies. The system described was coupled to a size-exclusion column to allow the monitoring of activity.

[0010] These prior art techniques which are of a sequential nature are not suitable for on-line coupling to liquid chromatographic or other fractionation systems, since these techniques do not allow the continuous monitoring of the fractionation effluent.

[0011] In the article of Irth et al. in the *Journal of Chromatography* **633** (1993) 65-72 and in the articles of Oosterkamp et al. in *Anal. Chem.* **66** (1994) 4295-4301, and in the *Journal of Chromatography* **653** (1994) 55-61, a method for the on-line detection of digoxigenin and its metabolites is described. The on-line detection process comprises the direct injection of a sample containing digoxigenin and its metabolites, a liquid chromatographic (LC) fractional separation step, the mixing of the effluent of the LC column with fluorescein-labeled antibodies against digoxigenin, the removal of free labeled antibodies from the mixture via passage through a small column packed with an antigen-bound support, and detection of the strongly fluorescent biochemical complexes.

[0012] The digoxigenin system described is based on association reactions of antibodies and antigens eluting from the analytical column. By the use of fluorescein-labelled antibodies detection limits in the nanomolar range are obtained.

[0013] This prior art assay is a heterogeneous detection system. It requires a separation step between free and bound label. Examples for free/bound label separation techniques are restricted-access phases and hollow-fibre modules.

[0014] The biochemical reagent in the previously described prior art assay consists of fluorescein-labelled fragments of anti-digoxigenin antibodies which were immunopurified and are commercially available. The commercial availability of purified, labelled antibodies is however exceptional. In almost all cases, antibodies are only available in unlabelled state in crudely purified antiserum. Although labelling and purification schemes for antibodies (or other affinity proteins) are known to the person skilled in the art, it will be preferred to use antisera, which may be commercially available, without any pretreatment.

[0015] Hsieh et al. (Molecular Diversity, 2 (1996), 189-196) describe a screening method using chromatography coupled with mass spectrometry. Their method requires a dual run, first in the absence of a library of protein targets and then with the library. This means that background compounds, especially the affinity proteins causing severe background signals over a broad range in the mass spectrum, are always present in this assay type, which limits the efficiency, selectivity and flexibility of the method considerably.

[0016] Kaur et al. (Journal of Protein Chemistry 16, No 5 (1997)) describe a method in which active components in a combinatorial library are bound to a receptor, isolated by size exclusion chromatography and trapped on a column after which the desorption solvent is used to dissociate the affinity protein/ligand complex. Subsequent identification is done by mass spectrometry, but with a background of the affinity protein.

[0017] The methodology described by Hsieh and Kaur is thus hampered by severe background signals in the mass spectrometer of the present affinity protein.

[0018] One of the advantages of the present invention is that it provides the possibility to detect the unknown ligand in the absence of a protein background in the MS. Especially for low concentration this is a requirement for successful assays.

[0019] One aim of the invention is to provide an improvement on the known assays by providing a technique, which does not require the use of known ligands.

[0020] The above mentioned problems are overcome by the provision of an on-line detection method comprising the on-line coupling of the effluent of a fractionation step to a mass spectrometer, which method comprises the addition of a controlled amount of an affinity molecule, such as an affinity protein, to said effluent, whereby the affinity molecules bind analytes in the effluent, followed by a separation step using a restricted-access support, whereby the analyte-affinity molecule complex is permeated, followed by a suitable dissociation step to dissociate the analyte-affinity molecule complex, followed by a second separation step in which the dissociated analyte and affinity molecules are separated, followed by detection of the analyte using the mass spectrometer.

[0021] The use of the first restricted access column provides the permeation of the analyte-affinity molecule complex while the smaller molecules are retained.

[0022] The dissociation step results in the dissociation of the complex. By subsequently separating the affinity molecules from the effluent the analyte can be directed to the MS.

[0023] In particular, a suitable affinity molecule for analytes to be detected - which affinity molecule may e.g. be an affinity protein - is added to the effluent of a liquid chromatography or a capillary electrophoresis column to react with the analytes eluting from the column.

[0024] The second separation step is carried out either by using a restricted-access support or by using a

hollow fiber support.

[0025] When using a restricted-access support in the second separation step the affinity molecule is retained, followed by elution of the affinity molecule from the restricted-access support using a suitable carrier stream, and directing the eluted stream to the mass spectrometer. Prior to elution of the affinity molecule, a washing step can be carried out to remove the dissociation solvent and the receptor proteins from the restricted-access column. In order to prevent false positives, it is important that the conditions, with respect to concentrations of compounds such as solvents, of the second restricted-access column are substantially constant and substantially the same as in the first restricted-access column.

[0026] When using a hollow fiber support in the second separation step the ligand is permeated and the permeate is directed to the mass spectrometer.

[0027] Preferably the dissociation step is a low pH shock, such as injection of a hydrochloric acid solution with a pH below 2.5, contacting with a high ionic strength solution, contacting with an organic solvent and/or contacting with a chaotropic reagent.

[0028] The fractionation step according to the present invention can be a liquid chromatography separation, a capillary electrophoresis step or a combinatorial chemistry system, which is optionally followed by a separation step which removes the high molecular weight background.

[0029] Preferred liquid chromatography fractionation steps comprise HPLC, reversed phase HPLC, capillary electrophoresis (CE), capillary electrochromatography (CEC), isoelectric focusing (IEF) or micellar electrokinetic chromatography (MEKC), all of which techniques are known to the person skilled in the art.

[0030] All known MS options and techniques can be used. Preferably the mass spectrometer is of the type chosen from the group consisting of electrospray ionization type, atmospheric pressure ionization type, quadrupole type, magnetic sector type, time-of-flight type, MS/MS, MS<sup>n</sup>, FTMS type, ion trap type and combinations thereof.

[0031] The present invention enables the detection of ligands for orphan receptors. Therefore the present invention is also directed to any such new ligand, found with the method set out in the above, as well as the use of these ligands as a ligand for orphan receptors.

[0032] The combination of the fractionation step with the MS detection step in accordance with the present invention greatly enhances the performance of both techniques. The combined techniques provide an analytical method, which is characterized by a high selectivity and a high sensitivity. Further, the problem associated with cross-reactivity does not occur when using the method of the present invention.

[0033] The method of the present invention uses bio-affinity molecules such as receptors, to detect any compounds showing high affinity for the ligand binding site

of said affinity molecule. The compounds to be detected may be biochemical compounds but are not in any way restricted thereto.

[0034] A system of high-throughput screening to be used in drug discovery may, for instance, comprise the following steps. Complex samples generated for instance by an upstream combinational chemistry system are prefractionated in fractions containing compounds of similar polarity using, e.g., a solid-phase extraction technique or electrophoretic sample handling principles. Each fraction may additionally be fractionated using, e.g., either analytical or preparative-scale liquid chromatographic fractionating columns. The compounds eluting from said LC column are on-line detected using a suitable affinity molecule detection technique. Where preparative-scale fractionating columns are applied, a postcolumn flow-split will be made. One of the two flow streams is subjected to detection using the affinity molecule detection technique; the other stream is directed to a fraction collector. Dependent on the signal obtained from the affinity molecule detector, fractions containing compounds causing a positive response will be collected while fractions causing a negative response will be discarded. This complete screening method can be automated using known valve-switching processes.

[0035] A suitable fractionation method to be used in the methods of the present invention comprises a liquid chromatography fractionation or a capillary electrophoresis step. Other fractionation techniques which are known to the person skilled in the art and which allow a relatively continuous output stream can, however, be used as well.

[0036] In a preferred embodiment, the liquid chromatography separation step is a reversed phase HPLC step.

[0037] It will be understood that according to the present invention as the effluent of a fractionation step is also to be understood the effluent of a flow injection, in which a mixture of different compounds is injected, optionally followed by a trapping step, for example using a restricted-access column. The separation step of the method according to the present invention then fulfills the fractionation, since the compounds of interest are fractionated by the separation step. The required selectivity can be obtained by operating the MS selectively tracing of the MS.

[0038] Since the flow injection technique provides a flexible and fast screening method, it is a preferred embodiment of the present invention.

[0039] All modes of MS-operation are possible in Flow Injection mode, however, due to the higher background, especially the very selective modes which comprise detecting ions of a selected single  $m/z$  trace or of selected multiple  $m/z$  traces are suitable.

[0040] Restricted-access supports are well-known to the person skilled in the art. In general, restricted-access supports, such as  $C_{18}$ -silica (for example Lichrospher™ C18 ADS from Merck) packed columns,

are capable of retaining smaller sized molecules while permeating the larger molecules. The type of restricted-access supports and the conditions under which they are operated, will depend on the type of assay. A person skilled in the art will be able to perform routine experiments to establish the optimal conditions and materials for each case.

[0041] Hollow fiber materials, for example made of polysulphon (obtainable from Amicon, Ireland) are well-known to the person skilled in the art. In general, hollow fiber modules are capable of permeating smaller sized molecules while retaining the larger molecules. Also the type of hollow fiber supports and the conditions under which they are operated, will depend on the type of assay. Again, the person skilled in the art will be able to perform routine experiments to establish the optimal conditions and materials for each case.

[0042] As already mentioned, the fractionation step can be a liquid chromatography separation, a capillary electrophoresis step or a combinatorial chemistry system. It has been found that it is very advantageous when this fractionation step is followed by a separation step, for example by using a hollow-fiber module, since this removes the high molecular weight background. Preferred liquid chromatography fractionation steps comprise HPLC, reversed phase HPLC, CE, CEC, IEF or MEKC, all of which techniques are known to the person skilled in the art.

[0043] All possible variations in MS techniques known in the art can be profited from according to the present invention. preferably, the MS is of the type chosen from the group consisting of electrospray ionization type, atmospheric pressure ionization type, quadrupole type, triple quadrupoles (QQQ) type, magnetic sector type, time-off-flight type, MS/MS,  $MS^n$ , FTMS type, ion trap type and combinations thereof.

[0044] For example, in scanning mode to trace compounds, low resolution MS with all possible instrumental designs of MS can be used, in particular quadrupole, magnetic sector, time-off-flight, FTMS and ion-trap. This generates typically molecular weight data with nominal mass accuracy.

[0045] When high resolution MS is applied, using all possible high resolution instrumental designs, in particular magnetic sector, time-off-flight, FTMS and ion trap, molecular weight data with high mass accuracy combined with the elemental composition of the compound can be obtained.

[0046] Other known MS techniques comprise tandem MS, such as MS/MS or  $MS^n$  (for example  $MS^3$ ). Application of these techniques enables the collection of structural information of the ligands, which is a preferred embodiment of the present invention. The data in scanning mode can be acquired in data-dependent mode which means that for each peak observed automatically the tandem MS measurement is performed. In addition to this the fragmentation induced in the tandem MS measurement can be coupled by more sophisticated

procedures, for example a broad band excitation, which also fragments expected fragments of less importance, such as the [protonated molecule-H<sub>2</sub>O]<sup>+</sup> peak in natural product screening.

**[0047]** The MS can also be operated in single/multiple ion monitoring mode, which can be used for highly selective detection.

**[0048]** Single/multiple ion monitoring mode can be used for highly selective detection. Using low resolution MS in single/multiple ion monitoring mode with all possible low resolution instrumental designs, in particular quadrupole, magnetic sector, time-off-flight, FTMS and ion trap, selective detection based on monitoring of previously determined m/z trace can be obtained. When high resolution MS is applied all possible instrumental designs of MS can be used, especially quadrupole, magnetic sensor, time-off-flight, FTMS and ion-trap. This enables typically very selective detection based on monitoring in high resolution previously determined m/z trace. When tandem MS is applied, all possible instrumental designs, in particular ion trap, quadrupole-time-off-flight (Q-TOF), FTMS and combinations of sector instruments with quadrupoles and ion traps, can be used. This enables the very selective detection based on monitoring a resulting peak or peaks in MS/MS and/or MS<sup>n</sup> measurements.

**[0049]** The assays according to the present invention can be based on continuous monitoring of the HPLC effluent after a separation step, for instance with a hollow fiber step, to remove the high molecular weight background. The MS is typically set to operate in one of single/multiple ions mode options.

**[0050]** With the assays according to the present invention the tracing of (bio-) active compounds in solutions of complex nature, for instance biological fluids or extracts, natural product extracts, solutions or extracts from biotechnological processes, resulting from chemical experiments, such as combinatorial technologies and processes and the like can be performed with a higher efficiency, selectivity and flexibility. Moreover, the present invention provides the possibility to limit the disturbance of background compounds.

**[0051]** The present invention further enables the identification of compounds based on conventional mass spectra, high resolution data or MS<sup>n</sup> based spectra.

**[0052]** With the method of the present invention it is also possible to perform library searching, based on a variety of mass spectrometric experiments enabling the screening of large series of samples and classifying these in classes based on similar active compounds, without the need for full identification of the compounds.

**[0053]** The assay method of the present invention can be applied in miniaturized formats of assay-based systems, for instance with chip technology based screening systems.

**[0054]** In the methods of the present invention all molecules capable of interaction with or binding to other molecules can be used as affinity molecule.

**[0055]** On-line coupling as used in the methods of the present invention requires fast reaction times in order to minimize extra-column band broadening. This means that affinity molecule-ligand interactions having reaction times in the order of minutes rather than hours should be considered. The suitable binding conditions under which the affinity molecules bind to the analyte comprise a contact time, which is in the same order of magnitude. Suitable binding conditions are conditions that provide optimal binding between the affinity molecules and the analyte. It will be understood that the precise conditions, such as temperature, residence time, chemical composition, will depend strongly on the type of assay and the proteins used therein.

**[0056]** The present invention will be described in further detail while referring to the following examples.

### Example 1 Restricted-access separation

**[0057]** Phosphate buffered saline (PBS) carrier phase solution (PBS), affinity protein solution, dissociation solution and desorption solution were pumped by Knauer ministar pumps. The flushing solution was added using a Minipuls 3 peristaltic pump. Two Rheodyne sixport valves were used to enable flushing and column switching. A third rheodyne injection valve was used to perform the loop injections (5 µl injection volume). The carrier phase consisted of PBS buffer (pH 7.0) which was pumped at a flow rate of 40 µl/min. A solution of anti digoxigenin Fab fragments (5 µM) in PBS was added to carrier phase at 80 µl/min. During control experiments the anti digoxin Fab fragments were switched with either PBS or PBS solutions containing 2 µM of goat antibiotin. The dissociation solution consisted of 1 M hydrochloric acid which was added at 260 µl/min. The desorption buffer consisted of 80% MeOH and 2mM ammonium acetate. The rinsing solution was composed of 5% isopropanol/2 mM ammonium acetate and was used to flush the second RA-column for one minute after trapping at 300 µl/min. The dimensions of the first and second restricted access columns equaled 4 x 15mm and 2 x 10 mm respectively. The RA-material used was Li-chrospher C18 ADS (Merck). Reactions were carried out in knitted open tubular Teflon reaction coils (0.5 mm ID), using reaction times of 30 s and 30 s for the first and second biochemical reaction respectively. Detection was carried out by a Finnigan LCQ mass spectrometer.

**[0058]** After injection in the carrier phase, digoxin was allowed to react for 30 s with anti digoxigenin Fab fragments. After reaction Fab bound digoxin was able to pass the first restricted access column, whereas excess unbound molecules were trapped onto the first RA-column. After passing the first RA-column the affinity complex was dissociated by a pH shock using a solution of hydrochloric acid. Again the mixture was allowed to react for 30 s after which the ligand could be trapped onto the second RA-column. Three minutes after injection,

the second RA-column was flushed for one minute, after which the second RA-column was switched into the desorption solution. Digoxin was desorbed and subsequently detected by the MS.

[0059] Figure 1 presents a typical MS spectrum of digoxin. Apart from the protonated molecule (780.8), an ammonium cluster can be observed at 797.9, whereas the m/z ratios of 762.8 and 651.1 represent the loss of respectively water and a sugar group. Wideband activation of the protonated molecule (780.8 +/- 18) at 45% collision energy fragmentates both the parent mass as well as the ammonium complex to 651.1, which can be observed as the base peak in MS/MS spectra under the conditions used (figure 2). During assay operation digoxin was detected using both MS and MS/MS (m/z 797.9 and 651.1 respectively).

[0060] Injecting 5, 50 and 500 ng digoxin in the assay system with and without anti digoxigenin Fab fragments showed that digoxin was completely trapped when the affinity protein was absent. However addition of anti digoxigenin Fab fragments showed increasing response of digoxin (measured at m/z 651.1), demonstrating the ability of the system to detect digoxin using the MS-based orphan assay.

[0061] In order to check the robustness of the ms-bcd system, the anti-digoxigenin Fab fragment solution was replaced by goat anti-biotin. Injection of 0, 5, 50 and 500 ng digoxin showed that the signal obtained for 500 ng is still comparable to the blank, meaning that the excess of digoxin is adequately trapped on the first restricted access column.

## Example 2 Hollow fiber separation

[0062] The set-up consisted of two Pharmacia (Uppsala, Sweden) P3500 pumps which delivered the mobile phase and affinity protein solution, one Gilson 306 LC-pump which delivered an acidic solution and two Gilson Minipuls 3 peristaltic pumps which respectively created backpressure at the retentate outlet and added phosphate buffer to the permeate phase. Mobile phase, affinity protein solution and acidic solution were pumped at a flow-rate of 0.2 ml/min. Acidic conditions in the permeate phase were neutralized by the addition of a phosphate buffer (100 µl/min). Injections were performed by an Aspec autosampler equipped with two Rheodyne six port valves (injection volume 20 µl). A Finnigan LQC mass spectrometer was used for detection.

[0063] 0.15 ml/min PBS carrier (136.9 mmol/l NaCl, 2.7 mmol/l KCl, 9.2 mmol/l Na<sub>2</sub>HPO<sub>4</sub> and 1.8 mmol/l K<sub>2</sub>PO<sub>4</sub>) comprising eighty randomly chosen drugs, each in the 1 µM range, was mixed with 10 nM goat anti-biotin (0.15 ml/min, obtained from Pierce, Rockford, Illinois, USA) by means of an inverted Y-piece. The mixture is allowed to react for 48 seconds in a reaction coil prior to the separation of unbound ligands and affinity complexes by a restricted-access column (C18 alkyl diol, Merck Lichrospher™ ADS, 15 x 4 mm). Subsequently

after passing the restricted access column the pH of the mobile phase is lowered to 2 by addition of hydrochloric acid (0.3 ml/min), causing the affinity complex to dissociate in a reaction coil within 10 s reaction time. On-line filtration by means of a hollow fiber module (polysulfon, obtained from Amicon, 0.5 mm I.D., 30 kDa cutoff, Module I.D. was 0.85 mm) then separates affinity molecules from ligands. The ligands could all subsequently be detected with the MS.

## Claims

1. On-line detection method comprising the on-line coupling of the effluent of a fractionation step to a mass spectrometer, which method comprises the addition of a controlled amount of an affinity molecule to said effluent, whereby the affinity molecules bind analytes in the effluent, followed by a separation step using a restricted-access support, whereby the analyte-affinity molecule complex is permeated, followed by a suitable dissociation step to dissociate the analyte-affinity molecule complex, followed by a second separation step in which the dissociated analyte and affinity molecules are separated, followed by detection of the analyte using the mass spectrometer.
2. On-line detection method according to claim 1, in which the second separation step is carried out using a restricted-access support, in which the affinity molecule is retained, followed by elution of the analyte from the restricted-access support using a suitable carrier stream, and directing the eluted stream to the mass spectrometer.
3. On-line detection method according to claim 1, in which the second separation step is carried out using a hollow fiber support, whereby the analyte is permeated and the permeate is directed to the mass spectrometer.
4. On-line detection method according to any of the preceding claims, in which the dissociation step is a low pH shock, contacting with a high ionic strength solution, contacting with an organic solvent and/or contacting with a chaotropic reagent.
5. Method according to any of the preceding claims in which the fractionation step is a liquid chromatography separation, a capillary electrophoresis step or a combinatorial chemistry system, which is optionally followed by a separation step which removes the high molecular weight background.
6. Method according to claim 5, in which the liquid chromatography separation step is a HPLC, a reversed phase HPLC, a CE, a CEC, a IEF or a MEKC

step.

7. Method according to any one of the preceding claims, in which the mass spectrometer is of the type chosen from the group consisting of electro-spray ionization type, atmospheric pressure ionization type, quadrupole type, magnetic sector type, time-off-flight type, MS/MS, MS<sup>n</sup>, FTMS type, ion trap type and combinations thereof. 5
8. Method according to any one of the preceding claims, in which the mass spectrometer is set to detect ions of a selected single m/z trace, selected multiple m/z traces, in scanning mode or any sequential mode. 10 15
9. The method of any of the preceding claims wherein the affinity molecule is an affinity protein.
10. The method of any of the preceding claims wherein the affinity molecule is an orphan receptor. 20
11. Compound detected by the method of any one of the preceding claims. 25
12. The use of a compound of the preceding claim as a ligand for affinity molecules. 30

35

40

45

50

55

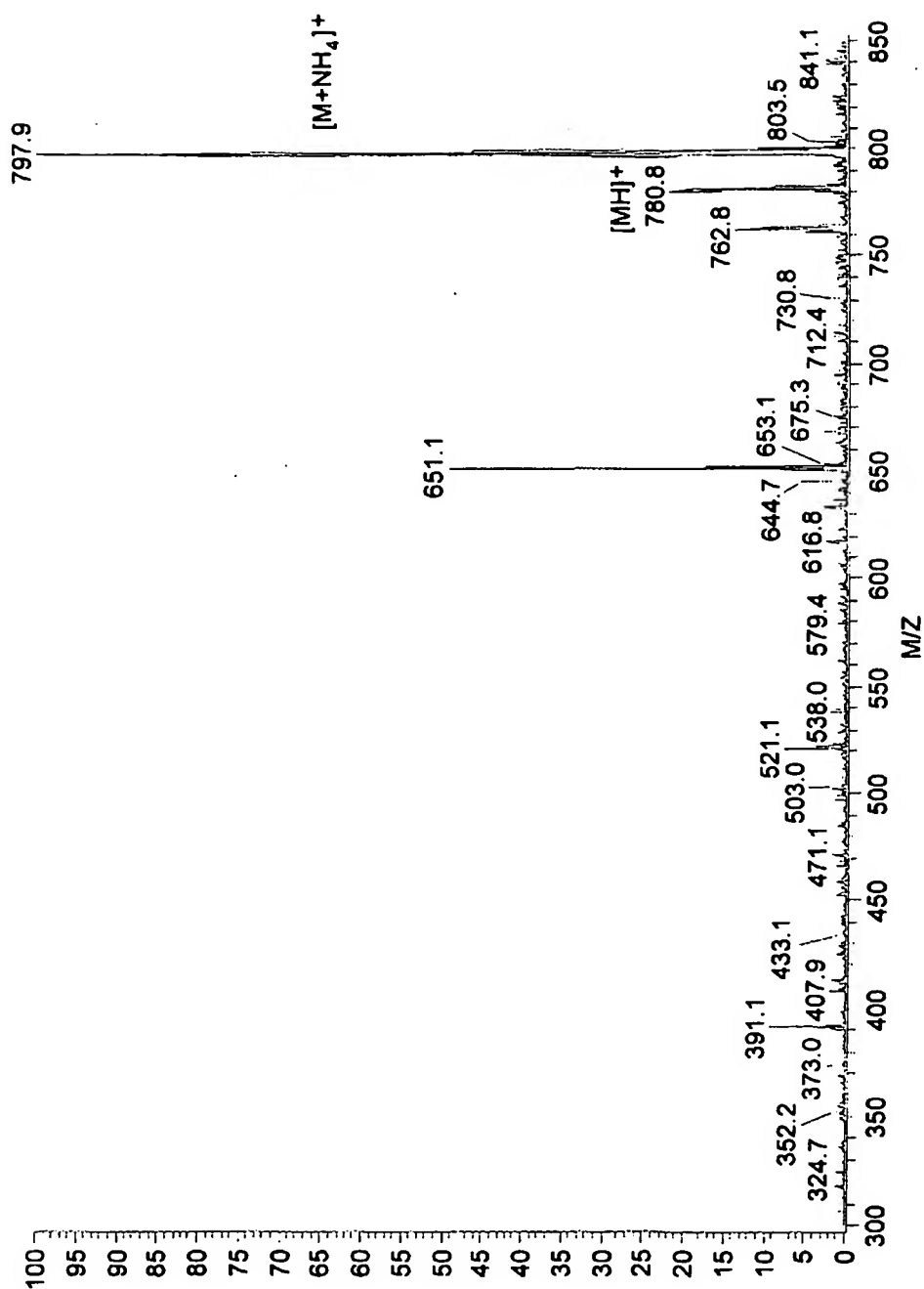


FIG. 1



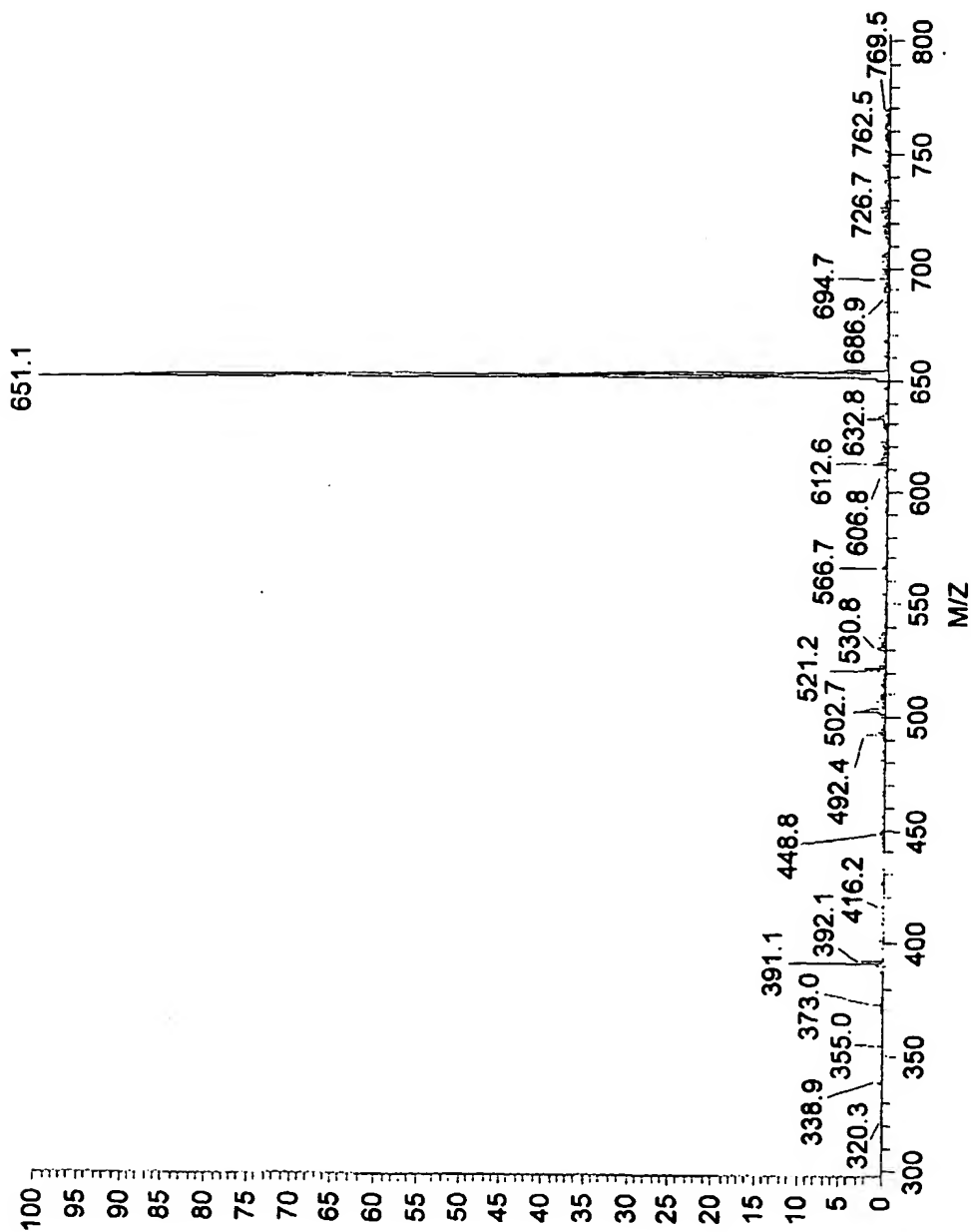


FIG. 2



European Patent  
Office

# PARTIAL EUROPEAN SEARCH REPORT

Application Number

which under Rule 45 of the European Patent Convention shall be considered, for the purposes of subsequent proceedings, as the European search report

EP 99 20 1301

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.7)
X	NEDVED, MICHAEL L.; HABIBI-GOUDARZI, SOHRAB; GANEM, BRUCE; HENION, : "Characterization of Benzodiazepine "Combinatorial" Chemical Libraries by Online Immunoaffinity Extraction, Coupled Column HPLC-Ion Spray Mass Spectrometry-Tandem Mass Spectrometry " ANALYTICAL CHEMISTRY, vol. 68, no. 23, 1 December 1996 (1996-12-01), pages 4226-4236, XP002117165 * page 4228 - page 4231 *	1,3-9	G01N33/566 G01N30/46
X	SIEGEL, MARSHALL M.; TABEI, KEIKO; BEBERNITZ, GERALDINE A.; BAUM, ELLEN Z.: "Rapid methods for screening low molecular mass compounds non-covalently bound to proteins using size exclusion and mass spectrometry applied to inhibitors of human cytomegalovirus protease" JOURNAL OF MASS SPECTROMETRY, vol. 33, no. 3, 1998, pages 264-273, XP002117166 * the whole document *	1,3-9	TECHNICAL FIELDS SEARCHED (Int.Cl.7)  G01N C12Q
-/--			
<b>INCOMPLETE SEARCH</b>			
<p>The Search Division considers that the present application, or one or more of its claims, does/do not comply with the EPC to such an extent that a meaningful search into the state of the art cannot be carried out, or can only be carried out partially, for these claims.</p> <p>Claims searched completely :</p> <p>Claims searched incompletely :</p> <p>Claims not searched :</p> <p>Reason for the limitation of the search:</p> <p>see sheet C</p>			
Place of search <b>THE HAGUE</b>		Date of completion of the search <b>1 October 1999</b>	Examiner <b>Hart-Davis, J</b>
<b>CATEGORY OF CITED DOCUMENTS</b> X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document		T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons & : member of the same patent family, corresponding document	

EPO FORM 1503 03/82 (P4/C07)



European Patent  
Office

**INCOMPLETE SEARCH  
SHEET C**

Application Number  
EP 99 20 1301

Claim(s) not searched:  
11-12

Reason for the limitation of the search:

Present claims 11-12 relate to a compound defined by reference to a desirable characteristic or property, namely its ability to bind to an (unnamed) orphan receptor.

The claims cover all compounds having this characteristic or property, whereas the application provides support within the meaning of Article 84 EPC and/or disclosure within the meaning of Article 83 EPC for only a very limited number of such compounds. In the present case, the claims so lack support, and the application so lacks disclosure, that a meaningful search over the whole of the claimed scope is impossible. Independent of the above reasoning, the claims also lack clarity (Article 84 EPC). An attempt is made to define the compound by reference to a result to be achieved. Again, this lack of clarity in the present case is such as to render a meaningful search over the whole of the claimed scope impossible.

## PARTIAL EUROPEAN SEARCH REPORT

Application Number  
EP 99 20 1301

DOCUMENTS CONSIDERED TO BE RELEVANT			CLASSIFICATION OF THE APPLICATION (Int.Cl.7)
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	
X	DUNAYEVSKIY, YURIY M.; LAI, JAN-JI; QUINN, CHERYL; TALLEY, FRANK; VOUIROS, PAUL: "Mass spectrometric identification of ligands selected from combinatorial libraries using gel filtration" RAPID COMMUNICATIONS IN MASS SPECTROMETRY, vol. 11, no. 11, 1997, pages 1178-1184, XP002117167 * the whole document *	1,3-9	
X	BLOM, KARL F.; LARSEN, BARBARA S.; MCEWEN, CHARLES N. : "Determining Affinity-Selected Ligands and Estimating Binding Affinities by Online Size Exclusion Chromatography/Liquid Chromatography-Mass Spectrometry" JOURNAL OF COMBINATORIAL CHEMISTRY, vol. 1, no. 1, January 1999 (1999-01), pages 82-90, XP002117168 * page 82 - page 85 *	1,3-9	TECHNICAL FIELDS SEARCHED (Int.Cl.7)
X	HUYER, GREGORY; KELLY, JOHN; MOFFAT, JASON; ZAMBONI, ROBERT; JIA, ZONGCHAO; GRESSER, MICHAEL J.; RAMACHANDRAN, CHIDAMBARAM: "Affinity selection from peptide libraries to determine substrate specificity of protein tyrosine phosphatases" ANALYTICAL BIOCHEMISTRY, vol. 258, no. 1, 1998, pages 19-30, XP002117169 * page 21 *	1,3-9	



European Patent  
Office

## PARTIAL EUROPEAN SEARCH REPORT

Application Number  
EP 99 20 1301

DOCUMENTS CONSIDERED TO BE RELEVANT			CLASSIFICATION OF THE APPLICATION (Int.Cl.7)
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	
D, A	HSIEH YF, GORDON N, REGNIER F, AFEYAN N, MARTIN SA, VELLA GJ: "Multidimensional chromatography coupled with mass spectrometry for target-based screening" MOLECULAR DIVERSITY, vol. 2, no. 4, April 1997 (1997-04), pages 189-196, XP002117170 * figures 1-7 *	1-9	
D, A	KAUR S, MCGUIRE L, TANG D, DOLLINGER G, HUEBNER V: "Affinity selection and mass spectrometry-based strategies to identify lead compounds in combinatorial libraries" JOURNAL OF PROTEIN CHEMISTRY, vol. 16, no. 5, July 1997 (1997-07), pages 505-511, XP002117171 * the whole document *	1-9	TECHNICAL FIELDS SEARCHED (Int.Cl.7)
A	WO 98 56028 A (AROLE INC ;KYRANOS JAMES N (US); LI YU TYSR (US)) 10 December 1998 (1998-12-10) * examples 1,2 *	1-9	
A	WO 97 01755 A (PERSEPTIVE BIOSYSTEMS INC) 16 January 1997 (1997-01-16) * claim 1; example 1 *	1-9	
A	WO 97 43301 A (CIBA GEIGY AG ;GOELLER CHRISTINE (US); KELLY MICHELE ANN (US); THO) 20 November 1997 (1997-11-20) * claim 1 *	1-9	
A	WO 96 37777 A (NELSON RANDALL W ;WILLIAMS PETER (US); KRONE JENNIFER REEVE (US)) 28 November 1996 (1996-11-28) * claims 1-3; examples 1-8 *	1-9	

EPO FORM 1503 03.82 (P04C10)

**ANNEX TO THE EUROPEAN SEARCH REPORT.  
ON EUROPEAN PATENT APPLICATION NO.**

EP 99 20 1301

This annex lists the patent family members relating to the patent documents cited in the above-mentioned European search report.  
The members are as contained in the European Patent Office EDP file on  
The European Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

01-10-1999

Patent document cited in search report		Publication date	Patent family member(s)		Publication date
WO 9856028	A	10-12-1998	AU	7713198 A	21-12-1998
WO 9701755	A	16-01-1997	EP	0835446 A	15-04-1998
			JP	11509314 T	17-08-1999
WO 9743301	A	20-11-1997	AU	2889097 A	05-12-1997
WO 9637777	A	28-11-1996	NONE		

EPO FORM P0459

For more details about this annex : see Official Journal of the European Patent Office, No. 12/82

**This Page is Inserted by IFW Indexing and Scanning  
Operations and is not part of the Official Record**

**BEST AVAILABLE IMAGES**

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ BLACK BORDERS
- ☐ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
- ☒ FADED TEXT OR DRAWING
- ☒ BLURRED OR ILLEGIBLE TEXT OR DRAWING
- ☐ SKEWED/SLANTED IMAGES
- ☐ COLOR OR BLACK AND WHITE PHOTOGRAPHS
- ☐ GRAY SCALE DOCUMENTS
- ☐ LINES OR MARKS ON ORIGINAL DOCUMENT
- ☒ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY
- ☐ OTHER: \_\_\_\_\_

**IMAGES ARE BEST AVAILABLE COPY.**

**As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.**